



KAPNOL MEDICAL
INSTRUMENTS

THE PODOLOGY CATALOG

Web: www.kapnol.com

Email: info@kapnol.com

Phone[#] (+92) 333 874 1182



*Photo of detail:
Fastening point of an
KAPNOL corner clipper*

A look at the manufacturing process

Depending on the model, an instrument needs 50 or more working steps, and approximately two thirds are carried out by hand. The joints are forged using resistance of up to 100 tons. This is followed by punching out, heat treatment, sawing and milling processes. Hardening at over 1,000°C sets the required strength of the instrument – for example, tweezers must be more flexible than scissors. At the same time, homogenous hardness provides an even distribution of the chrome components and corrosion resistance from the surface to the core.

Polishing evens out the surface and the corrosion-resistant chrome oxide protection layer (passive layer) is finally activated. Instruments generally have a non-reflective matt finish.

The manufacturing experts (tool mechanics specializing in instrument engineering with three and a half years of training) may handle a complex tool such a nail clipper approx. 150 to 180 times before it leaves the production process.

I'll have stainless steel!

With the correct care, high quality KAPNOL instruments are extremely hard-wearing.

- With slight mechanical superficial damage there will not be any corrosion.
- Correctly cared for, the surfaces will remain intact and beautiful for years.
- Formation of rust on cutting surfaces and rust beds at joints and friction surfaces can be avoided in the long term.

Important for hygiene: even when stainless steel instruments are placed in cleaning and disinfectant solutions several times a day, this does not cause "flaking" or other changes to the surface of the material.

And: in contrast to rust-forming metals and alloys, stainless steel can be sterilized in an autoclave. This means that the instruments can be safely sterilized and are quickly reusable.



Quality from the start

*Our instruments are highly valued by surgeons all over the world.
I know why this is the case and will ensure that things remain this way*

Golden rules for instrument care


- Only use suitable disinfectants and cleaning agents to disinfect and clean instruments. Please use no drill brushes for manual cleaning. We recommend a medium hard toothbrush.
- When placing the instruments in the ultrasound bath, always ensure that the joints are opened and the joint area of the instruments is cleaned.
- When using ultrasound, take the immersion times and dosage of the cleaning agent into account.
- Prior to sterilization, rinse the instruments with fully desalinated or distilled water and dry them.
- Regularly spray the joint area of instruments with a special paraffin oil to prevent it from setting in.
- Absolutely sort out instruments which have rust or other safety-relevant damage.

This avoids costly instrument repairs

KAPNOL quality is reliable and permanent.
Over 20 years of expertise and care of the instrument
mechanics contribute to this.

Barrel sprung pliers

Protect against excessive
stress on the podiatrist's
hand and elbow joints.

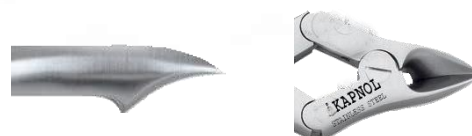

Nail clipper

For cutting the free edges of the nails.
Particularly stable, strong model.



Left side

Order No. MD

1 Nail clipper

15 cm, barrel sprung

Stainless steel, KAPNOL

HF 250 [309 200 000](#)



2 Nail clipper

12 cm, barrel sprung

Stainless steel, KAPNOL

G 495 [309 200 100](#)

4 HF 203 Nail clipper



12 cm, matted, cutting edge 11 mm

Stainless steel, KAPNOL

G 345 [309 200 200](#)

3 Nail clipper

Right side

4 Nail clipper

13 cm, applied,
curved cutting edge 22 mm

HF 203 [309 100 400](#)



5 HF 241 Nail clipper



3 Nail clipper

12 cm, applied,
curved cutting edge 16 mm

HF 241 [309 134 200](#)



6 HF 243 Nail clipper



6 Nail clipper

14 cm, applied,
curved cutting edge 19 mm

Stainless steel, KAPNOL

HF 243 [309 134 400](#)

Nail clipper

For cutting the free edges of the nails. Particularly stable, strong model.



HF 210 Nail clipper



HF 211 Nail clipper



HF 223 Nail clipper



Left side

Order No. MD

1 Nail clipper

(head cutter) 13.5 cm, box jointed, cantilevered cutting edge 20 mm

Stainless steel, KAPNOL

HF 210 309 111 200

2 Nail clipper

(head cutter) 13.5 cm, box jointed, cantilevered cutting edge 14 mm

Stainless steel, KAPNOL

HF 211 309 111 300

3 Nail clipper

13 cm, box jointed, straight cutting edge 20mm

Stainless steel, KAPNOL

HF 223 309 122 400

Right side ►

4 Corner clipper

11.5 cm, straight cutting edge 15 mm

Stainless

HF 475 309 187 800

5 Corner clipper

11.5 cm, extra slim, straight cutting edge 13mm

Stainless steel, KAPNOL

HF 476 309 187 900

6 Corner clipper

11.5 cm, straight cutting edge 15 mm

Stainless steel, KAPNOL

HF 477 309 188 000

7 Corner clipper

11.5 cm, straight cutting edge 16 mm

Stainless steel, KAPNOL

HF 481 309 186 700

Corner clipper (small)

For cutting the sides and the free distal edges of the nails. For removal of the ingrown part of the nail (Unguis incarnatus) in the nail fold or for the removal of hyperkeratosis and clavi (corns and hard skin).



4 HF 475 Corner clipper



5 HF 476 Corner clipper



7 HF 481 Corner clipper



Corner clipper (large)

For cutting the sides and the free distal edges of the nails and thin soft nail plates. For removal of the ingrown part of the nail (Unguis in- carnatus) in the nail fold or for the removal of hard skin and calluses



1 HF 478 Corner clipper



2 HF 479 Corner clipper



3 HF 480 Corner clipper



Left side

Order No. MD

13 cm, straight cutting edge 17 mm

Stainless steel, KAPNOL

HF 478 [309 188 100](#)

2 Corner clipper

13 cm, straight cutting edge 17 mm

Stainless steel, KAPNOL

HF 479 [309 186 000](#)

3 Corner clipper

15 cm, straight cutting edge 18 mm

Stainless steel, KAPNOL

HF 480 [309 186 100](#)



4 HF 483 Corner clipper



Right side

Corner clipper

13 cm, straight cutting edge 17 mm

Stainless steel, KAPNOL

HF 483 [309 186 200](#)

5 Corner clipper

12 cm, cutting edge 16 mm

Stainless steel, KAPNOL

HF 482 [309 187 000](#)

6 Corner clipper

11.5 cm, straight cutting edge 13.5 mm

Stainless steel, KAPNOL

G 490 [309 185 100](#)



5 HF 482 Corner clipper



6 G 490 Corner clipper





1 G 491 Corner clipper



2 HF 470 Corner clipper



3 HF 460 Skin clipper



4 HF 461 Skin clipper



Skin clipper

For the removal of all kinds of clavi from the feet, such as clavi between the toes, hyperkeratosis and protruding edges of corns.

Left side

Order No. MD

10 cm, straight cutting edge 11 mm

Stainless steel, KAPNOL

G 491 [309 185 000](#)

Corner clipper

10 cm, straight cutting edge 11 mm

Stainless steel, KAPNOL

HF 470 [309 188 200](#)

3 Skin clipper

10 cm, cantilevered edge 7 mm

Stainless steel, KAPNOL

HF 460 [309 166 200](#)

11.5 cm, cantilevered edge 7 mm

Stainless steel, KAPNOL

HF 461 [309 166 300](#)

Right side

11,5 cm, cantilevered edge 11 mm

Stainless steel, KAPNOL

HF 465 [309 166 700](#)

Nail skin clipper

10 cm, cutting edge 10 mm

Stainless steel, KAPNOL

HF 448 [309 188 300](#)

Nail splitter, 13.5 cm

Stainless steel, K A P N O L

HF 266 [309 193 300](#)



5 HF 465 Skin clipper




6 HF 448 Nail skin clipper




7 HF 266 Nail splitter




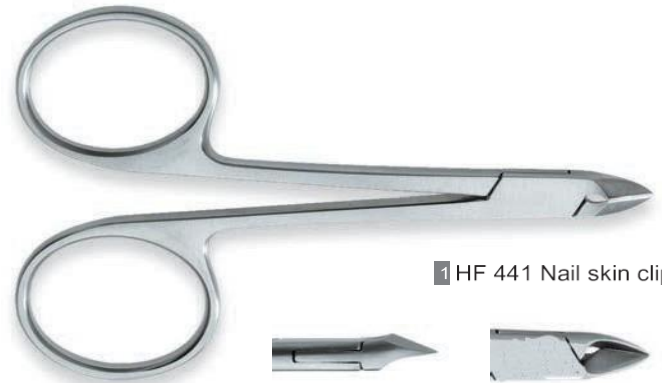


1 HF 441 Nail skin clipper

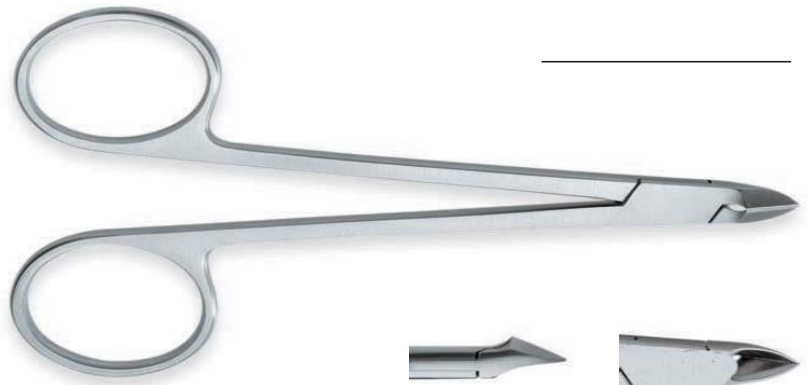

2 HF 443 Nail skin clipper

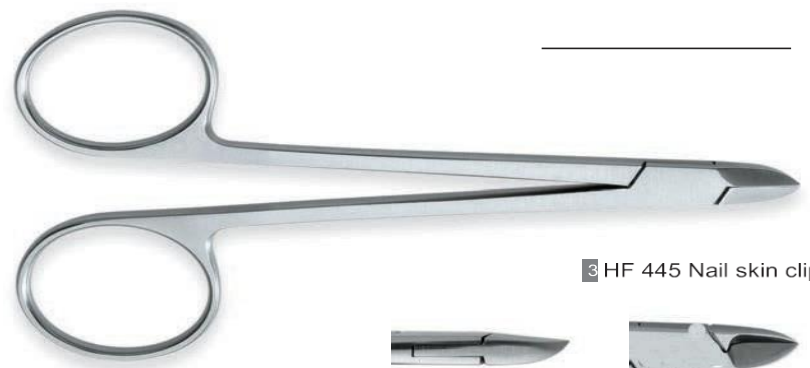

3 HF 445 Nail skin clipper

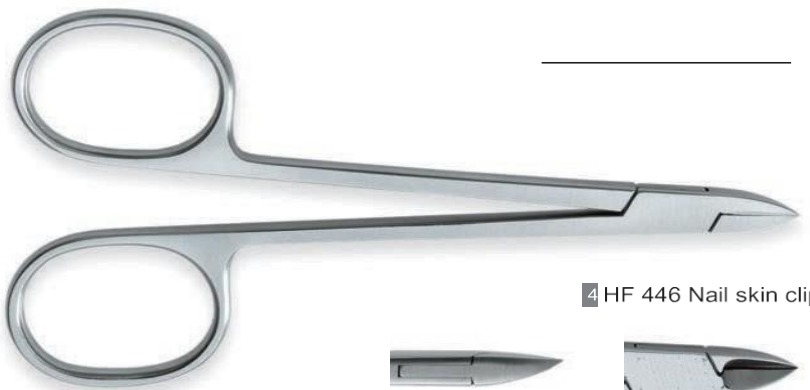

4 HF 446 Nail skin clipper

◀ Left side

Order No. MD

1 Nail skin clipper

8.5 cm, curved cutting edge 7 mm

Stainless steel, KAPNOL

HF 441 [309 154 200](#)

2 Nail skin clipper

10.5 cm, curved cutting edge 8 mm

Stainless steel,

3 Nail skin clipper

10.5 cm, bevelled cutting edge 8 mm

Stainless steel, KAPNOL

HF 445 [309 154 600](#)

4 Nail skin clipper

10.5 cm, straight cutting edge 8 mm

Stainless steel, KAPNOL


HF 446 [309 154 700](#)



5 BC 324 Dressing scissor (Illustrated smaller than original size)



6 HF 30 Nail scissor



7 HF 31 Nail scissor



8 HF 409 Nail skin scissor

Right side

5 Dressing scissor

14.5 cm, straight cutting edge 48 mm blunt-tipped

Stainless steel, KAPNOL



[309 192 400](#)

6 Nail scissor

9.5 cm, curved cutting edge Stainless steel,

KAPNOL

Stainless steel, KAPNOL

HF 30 [309 143 000](#)

6 Nail scissor

10.5 cm, curved cutting edge

Stainless steel, K A P N O L

Stainless steel, KAPNOL

HF 31 [309 143 100](#)

8 Nail skin scissor

9 cm, curved cutting edge 16 mm

Stainless steel, KAPNOL

Stainless steel, KAPNOL

HF 409 [309 141 000](#)



KAPNOL MEDICAL INSTRUMENTS

Tweezers Removal of hard skin, nail parts and foreign bodies.
Applying treatment materials in the nail fold.
Applying and removing wound dressings.

Order No. MD



OC 069R Anatomical tweezers

Anatomical tweezers, 12 cm
Stainless steel, KAPNOL
OC 069R [310 112 500](#)



BD 223R Anatomical tweezers

Anatomical tweezers, 15 cm
Stainless steel, KAPNOL
BD 223R [310 112 700](#)



BD 303 Splinter tweezers

Splinter tweezers, 11.5 cm
straight
Stainless steel, KAPNOL
BD 303 [310 106 600](#)



BD 023 Anatomical tweezers

Anatomical tweezers, 11.5 cm
moderately wide
Stainless steel, KAPNOL
BD 023 [310 112 300](#)



OC 220 Anatomical tweezers

Anatomical tweezers, 12 cm
Stainless steel
OC 220 [310 112 600](#)



KAPNOL MEDICAL INSTRUMENTS

Order No. MD

Splinter tweezers, 8 cm

Stainless steel, KAPNOL

BD 963 [310 106 300](#)



BD 963 Splinter tweezers

Splinter tweezers, 9.5 cm

Stainless steel, KAPNOL

BD 964 [310 106 400](#)



BD 964 Splinter tweezers

Splinter tweezers, 11.5 cm

Stainless steel, KAPNOL

BD 965 [310 106 500](#)



BD 965 Splinter tweezers

Splinter tweezers, 11.5 cm

Stainless steel, KAPNOL

BD 312 [310 101 200](#)



BD 312 Splinter tweezers

Eyelash tweezers, 8.5 cm

Stainless steel, KAPNOL

HH 225 [310 120 000](#)



HH 225 Eyelash tweezers

Nail Clippers

Stainless steel instruments in the best craftsmanship quality

Stainless steel instruments

These stainless steel instruments are produced based on many years of experience according to the best Surgical tradition in Pakistan. These high quality instruments are produced in many work steps on modern machines and the best manual work. Many years of tradition and experienced, well trained employees as well as using high quality steels and materials lay the foundation for fulfilling all medicinal product guideline requirements and being allowed to bear the CE label without restrictions. Constant material tests and strict multiple controls as part of quality assurance guarantee consistent quality. With proper handling and care, these tested and certified instruments will please you for a long time.

Nail clipper

Trimming the open nail edges.
Especially sturdy, strong models.



1 234 Nail clipper



3 243 Nail clipper



Left side

Order No. MD


14 cm, applied,
cantilevered cutting edge 22 mm
Stainless steel
234 [309 210 500](#)

13 cm, applied,
cantilevered cutting edge 22 mm
Stainless steel
203 [309 210 100](#)

13.5 cm, applied,
cantilevered cutting edge 19 mm
Stainless steel
243 [309 210 700](#)

Nail clipper

Trimming the open nail edges.
Especially sturdy, strong models.

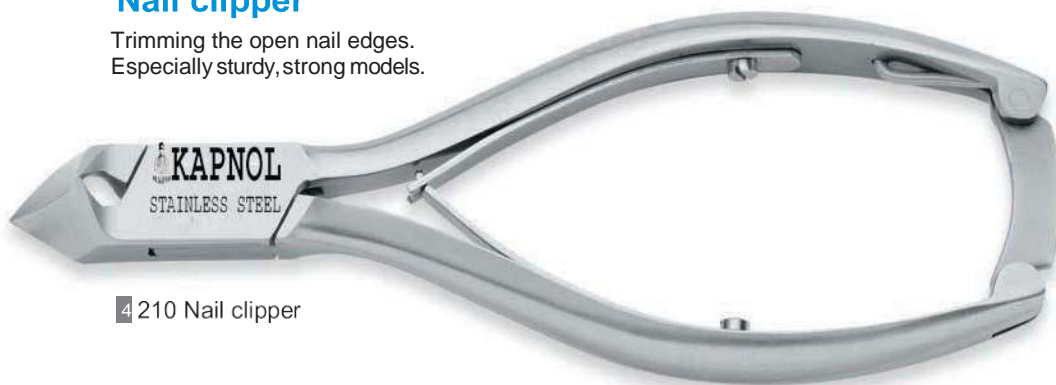

4210 Nail clipper



495 Nail clipper



Right side

(head cutters)

13.5 cm, box jointed,
cantilevered cutting edge 20 mm
Stainless steel
210 [309 210 200](#)

12 cm, barrel sprung,
Stainless steel
495 [309 250 300](#)



6211 Nail clipper



(head cutters)

13.5 cm, box jointed,
cantilevered cutting edge 14 mm
Stainless steel
211 [309 210 300](#)

Clippers

Nail clipper

Trimming the open nail edges.
Especially sturdy, strong models.



1 345 Nail clipper



Corner clipper

Cutting the lateral, distally located
open nail edges and thinner soft
nail plates. Removing ingrown
nail segments (unguis incarnatus)
in the nail edge and removing
callus and weals.



2 475 Corner clipper



3 476 Corner clipper



4 481 Corner clipper



KAPNOL MEDICAL INSTRUMENTS

Left side

Order No. MD

12 cm, applied,
cantilevered cutting edge 11 mm
Stainless steel
345

309 213 400



5478 Corner clipper

Corner clipper
11.5 cm, straight cutting edge 15 mm
Stainless steel
475

309 212 300



Corner clipper
extra slim, 11.5 cm,
straight cutting edge 14.5 mm
Stainless steel
476

309 212 400



6479 Corner clipper

Corner clipper
11.5 cm, straight cutting edge 16 mm
Stainless steel
481

309 212 100



Right side

Corner clipper
13 cm, straight cutting edge 17 mm
Stainless steel
478

309 212 600



7480 Corner clipper
(Illustrated smaller than original size)

Corner clipper
13 cm, straight cutting edge 17 mm
Stainless steel
479

309 211 700



Corner clipper
15 cm, straight cutting edge 18 mm
Stainless steel
480

309 211 900

Clippers

Corner clipper (small)

Cutting the lateral, distally located open nail edges. Removing ingrown nail segments (unguis incarnatus) in the nail edge and removing hyperkeratoses and clavi (corns and callus).



1 483 Corner



2 486 Corner clipper



3 490 Corner clipper



4 491 Corner clipper



Left side

Order No. MD

Corner clipper

13 cm, straight cutting edge 17 mm

Stainless steel

483 309 212 000



5 492 Corner clipper

Corner clipper

11.5 cm, straight cutting edge 24.5 mm

Stainless steel

486 309 213 600



Corner clipper

11.5 cm, straight cutting edge 13.5 mm

Stainless steel

490 309 213 200



5 470 Corner clipper

Corner clipper

10 cm, straight cutting edge 11 mm

Stainless steel

491 309 213 100



Right side

Corner clipper

11.5 cm, straight cutting edge 13.5 mm

Stainless steel

492 309 213 700



8 2510 Nail and Corner clipper

Corner clipper

10 cm, straight cutting edge 11 mm

Stainless steel

470 309 212 700



Corner clipper

11.5 cm, straight cutting edge 15mm

Stainless steel

477 309 212 500



7 477 Corner clipper

Nail and Corner clipper

13 cm, cantilevered cutting edge 14

Stainless steel

2510 309 251 000



Clippers

Skin clipper

For the removal of all kinds of clavi from the feet, such as clavi between the toes, hyperkeratosis and protruding edges of corns.



1 460 Skin clipper



2 461 Skin clipper



3 465 Skin clipper



4 448 Nail skin



Left side

Order No. MD

10 cm, cantilevered cutting edge 7 mm

Stainless steel

460 [309 211 400](#)

11.5 cm, cantilevered cutting edge 11 mm

Stainless steel

461 [309 211 500](#)

11.5 cm, cantilevered cutting edge 11 mm

Stainless steel

465 [309 211 600](#)

Nail skin clipper

10 cm, cutting edge 10 mm

Stainless steel

448 [309 212 800](#)



5 441 Nail skin clipper



Right side

Nail skin clipper

8.5 cm, cantilevered cutting edge 7 mm

Stainless steel

441 [309 211 000](#)



6 443 Nail skin clipper



Nail skin clipper

10 cm, cantilevered cutting edge 8 mm

Stainless steel

443 [309 211 100](#)



7 445 Nail skin clipper



Nail skin clipper

10.5 cm, bevelled cutting edge 8 mm

Stainless steel

445 [309 211 200](#)

Clippers



446 Nail skin clipper



2 409 Nail skin scissor



3 30 Nail scissor





Nail Splitter and Tweezers Skin Clipper

Left side

Order No. MD

Nail skin clipper

10 cm, straight cutting edge 8 mm

Stainless steel
446 309 211 300

Nail skin scissor

9 cm, curved cutting edge 16 mm

Stainless steel
409 309 210 800

Nail scissor

9.5 cm, curved cutting edge

Stainless steel
30 309 210 900



Nail splitter



Eye Brow Tweezers R

right side

Nail splitter

With replacement blade
Stainless steel

312 203 300

Replacement blades for nail splitter

Without illustration 312 203 400

Tweezer skin clipper

Stainless steel



Diamond bar kit

Eye brow tweezer R


9 cm, stainless

310 115 000

Diamond bar kit

310 116 100

Scalpels and Blades

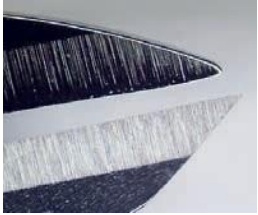

KAPNOL scalpel blades
with improved point



BB063R - for blades BB 10-15



BB 73 - for blades BB 6-15



BB 84 - for blades BB 18-23



BB 6



BB 10



BB 11



BB 13



BB 18



BB 19



BB 21



BB 20



BB 22



BB 23

Scalpel grip KAPNOL

Description	Order No.	Amount	MD
BB063R	311 100 300	1 Item	
Round grip; Stainless steel, 14.5 cm, for blades BB 10-15			
Stainless steel, 12.5 cm, for blades BB 6-15	311 100 100	1 Item	
Stainless steel, 13.5 cm, for blades BB 18-23	311 100 200	1 Item	

KAPNOL sterile scalpel blades

Description/Order No.	Amount	MD
BB 6 311 100 600	10 Items 100 Items	
BB 10 311 101 000	10 Items 100 Items	
BB 11 311 101 100	10 Items 100 Items	
BB 13 311 101 300	10 Items 100 Items	
BB 15 311 101 500	10 Items 100 Items	
BB 18 311 101 800	10 Items 100 Items	

KAPNOL sterile scalpel blades

Description/Order No.	Amount	MD
BB 19 311 101 900	10 Items 100 Items	
BB 20 311 102 000	10 Items 100 Items	
BB 21 311 102 100	10 Items 100 Items	
BB 22 311 102 200	10 Items 100 Items	
BB 23 311 102 300	10 Items 100 Items	

In podiatry, specialist careful removal of hard skin and corns on areas of the feet that are especially subject to stress is part of the daily routine.

Use of the right scalpels, blades and handy blade holders makes treatment fast and pain free. Choose from the available models.

Replacement scalpels
for the removal of
corns
hard skin
calluses



Removal of corns with a scalpel.



Using light pressure, the scalpel is guided over tautened skin.



Removal of hard skin from the forefoot using a scalpel.



01001/A/1 - for blades 01001/00-6

Diabetic plier HF 99

Diabetic Plier

12 cm
Stainless steel KAPNOL
HF 99

Description: This plier is use for diabetic patients



Corner Clipper HF 100

Corner clipper

Size 13 cm,

Stainless steel KAPNOL
HF 100



Scalpels and Blades

KAPNOL flat hard skin rasp

Size 10 cm, Cutting edge 7 mm
Order No. HF 101


Nail skin clipper HF 101



Nail skin clipper

Size 10 cm, Cutting edge 9 mm
Size 12 cm, Cutting edge 10 mm
Stainless Steel KAPNOL
HF 102


Nail skin clipper HF 102

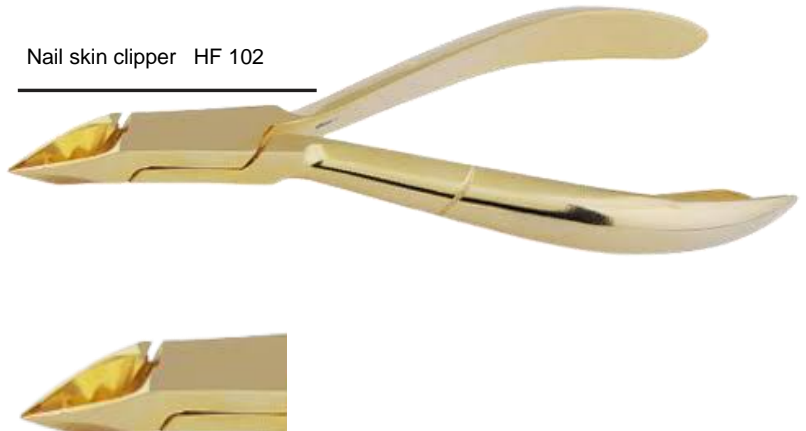

Nail clipper

Size 12 cm
Size 10 cm
Stainless Steel KAPNOL
HF 103

Nail Clipper HF 103

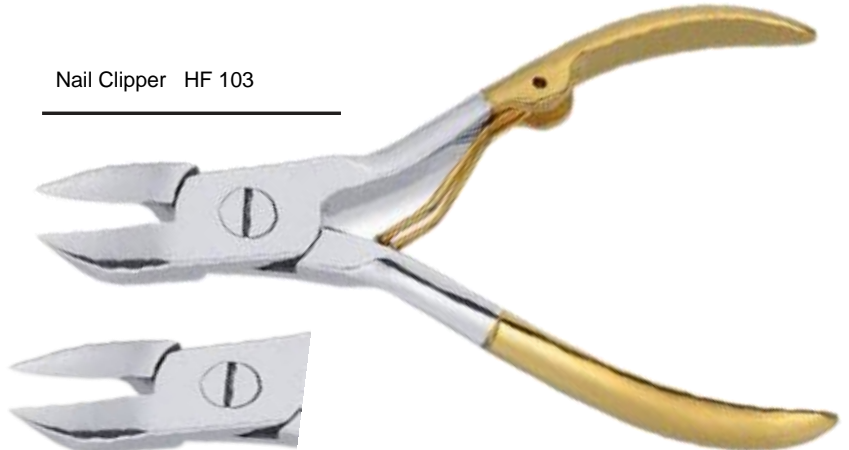

Nail skin clipper HF 104

Nail skin clipper

12 cm, cutting edge 8 mm
Stainless Steel KPNOL
HF 104



Nail skin clipper HF 105

Nail skin clipper

11.5Cm, cutting edge 9 mm
Stainless Steel KAPNOL
HF 105



Nail skin clipper

12 cm, Cutting edge 5 mm
Stainless Steel KAPNOL
HF 106




Nail skin clipper HF 106

Nail skin clipper

S 11 cm, Cutting edge 8 mm
Stainless Steel KAPNOL
HF 107

Nail skin clipper HF 107

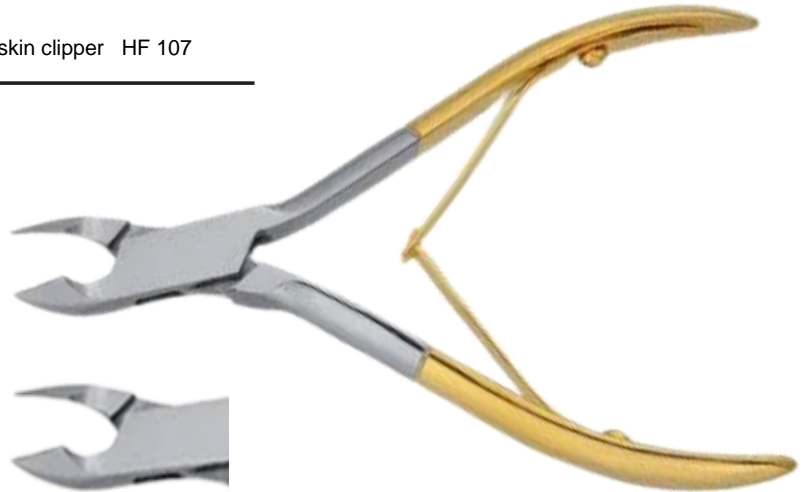

Stork Scissor HF 108

Stork Scissor


Size 10 cm
HF 108



Hollow blade holder

Size Standard
Order No. HF109

Hollow blade holder HF 109

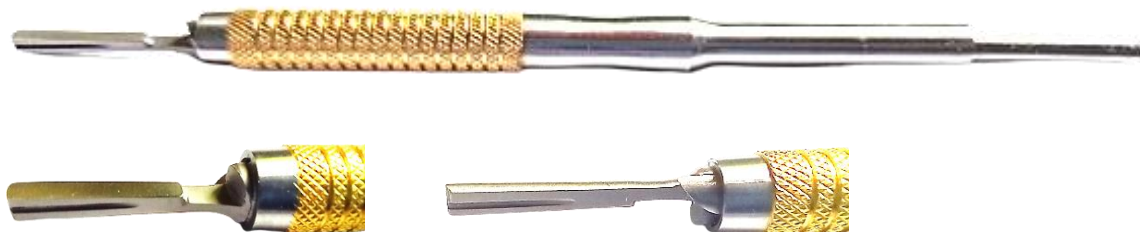

Hollow blade holder

Size Standard
Order No. HF 110

Hollow blade holder HF 110



KAPNOL

Size Standard
Order No. HF 111

Double instrument HF 111



Nail clipper

15 cm
Stainless steel KAPNOL
HF 112



Nail Clipper HF 112

Nail clipper

12 Cm , 14cm
Stainless steel KAPNOL
HF 113



Nail Clipper HF 113

Nail clipper

12 Cm , 14 Cm
Stainless steel KAPNOL
HF 114



Nail Clipper HF 114

Nail clipper

12 cm
Stainless steel KAPNOL
HF 115



Nail Clipper HF 115

Nail skin clipper

11.5 cm cutting edge 7 mm,
Stainless steel KAPNOL
HF 116



Nail skin Clipper HF 116

Nail skin clipper

11.5 cm cutting edge 6 mm
Stainless steel
HF 117



Nail Clipper HF 117

Nail skin clipper

10 cm cutting edge 7 mm
Stainless steel KAPNOL
HF 118



Nail Skin Clipper HF 118

Nail skin clipper

12 cm , cutting edge 7 mm
Stainless steel KAPNOL
HF 119



Nail skin Clipper HF 119

Nail skin clipper

12 cm cutting edge 9 mm
Stainless steel KAPNOL
HF 120



Nail skin clipper HF 120

Nail skin clipper

12 cm cutting edge 10 mm
Stainless steel KAPNOL
HF 121

Nail skin clipper HF 121



Nail skin clipper

10 cm cutting edge 9 mm
Stainless steel KAPNOL
HF 122

Nail skin clipper HF 122



Tip Cutter HF 123

Tip Cutter HF 123

Stainless steel KAPNOL
HF 123



Nail skin scissor

9 cm
Stainless steel KAPNOL
HF 124

Nail skin scissor HF 124



Nail skin scissor

9 cm
Stainless steel KAPNOL
HF 125

Nail skin scissor HF 125



Ingrown Nail Nipper

12 cm
Stainless steel KAPNOL
HF 126

Ingrown Nail Nipper HF 126






KAPNOL MEDICAL INSTRUMENTS

KAPNOL Forceps

Size 5 inch
Stainless Steel KAPNOL
Order No. HF 127

Forceps HF 127

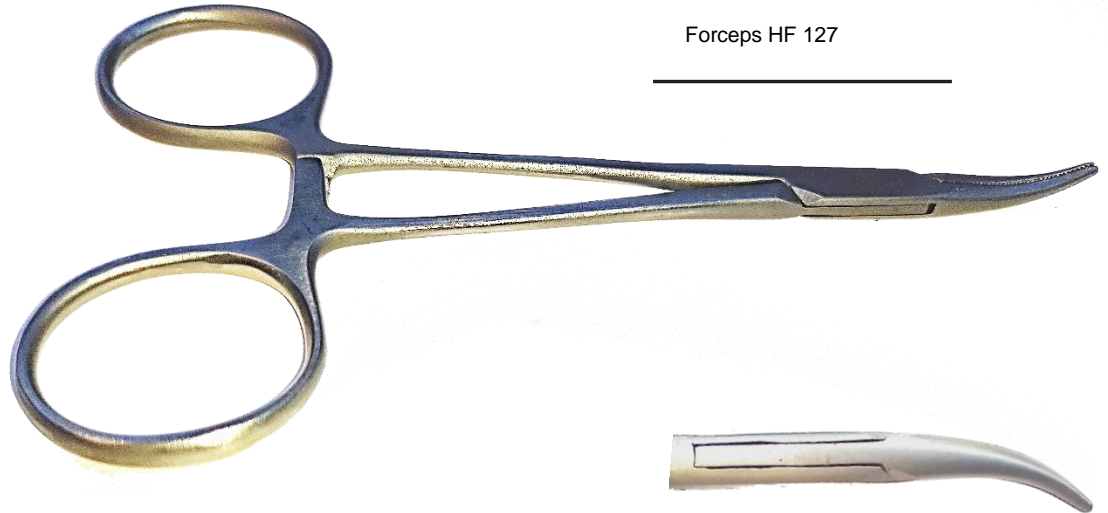

KAPNOL flat hard skin rasp


Size
Order No.




KAPNOL Fingers Separators

Size
Order No.

Fingers Separators HF 128

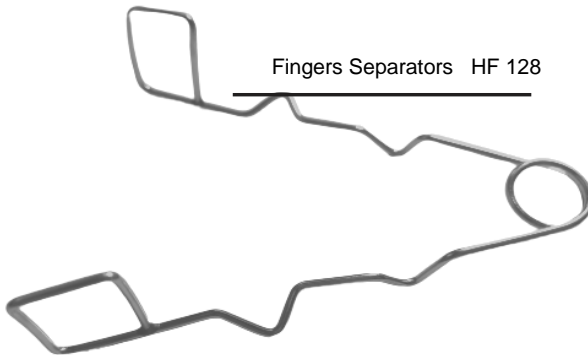

Fingers Separators HF 129

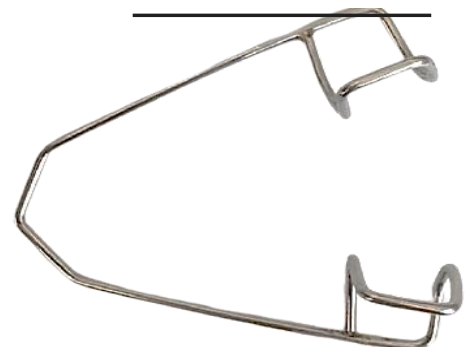



Hollow blade holder KEHR



Hollow blade holder CAPRON

Description

Hollow blade holder

Set, 3 Items

Blades CAPRON

Order No. MD

311 140 000

Blades PARAMOUNT

Order No. MD

311 160 000



Order No.

311 160 111



Order No.

311 160 112

Specialist Instruments and Files

Specialist instruments and files

For working on the nail,
under the nail and in the
nail fold.



Double instrument HH 106 R



Probes DC 64 R



Probes XG 056/WP 333



Double instrument DC 271 R



Probes DC 200 R



Probes 87035-160



Double spoons DC 241 R



Double spoons 87034-160



Hard skin hook HH 120 R

Left side

Order No.

Double instrument, 14.5 cm
Stainless steel, KAPNOL
HH 106 R [312 120 600](#)

Probes, 15 cm
Stainless steel, KAPNOL
DC 64 R [312 176 400](#)

Probes, 18 cm
Stainless steel, KAPNOL
XG 056 / WP 333

[312 173 000](#)

Double instrument, 17 cm
Stainless steel, KAPNOL
DC 271 R [312 120 700](#)

Probes, 17 cm
Double ended, Stainless
steel, KAPNOL
DC 200 R [312 172 000](#)

Probes, 16 cm
Stainless steel
87035-160 [312 210 100](#)

Double spoons, 17 cm
Stainless steel, KAPNOL
DC 241 R [312 143 800](#)

Double spoons, 16 cm
Stainless steel
87034-160 [312 210 200](#)

Hard skin hook, 15 cm
blunt, Stainless steel, KAPNOL
HH 120 R [312 132 000](#)



Spoon HH 118 R



Corner file HH 125 R



Corner file O 33 R



Corner file G 540 R



Corn chisel HH 60/1 R



Hoof pick

Right Side

Spoon, 12.5 cm
Stainless steel, KAPNOL
HH 118 R [312 143 700](#)

Corner file, 13.5 cm
Stainless steel, KAPNOL
HH 125 R [312 152 500](#)

Corner file, 13 cm
double ended, Stainless steel
O 33 R [312 153 300](#)

Corner file, 17.5 cm
double ended, Stainless steel,
KAPNOL
G 540 R [312 153 400](#)

Corn Chisel, 14 cm
1.8 mm, Stainless steel
HH 60/1 R [312 116 000](#)

Hoof pick, 11 cm

[318 010 000](#)

Special instruments, Files and Rasps

Rasps and files

Removal of hard skin and calluses.

Hard skin file stainless steel



Double hard skin file SN 32





KAPNOL MEDICAL INSTRUMENTS

Left side

Order No.

KAPNOL ceramic hard skin rasp

Wooden hard skin files

Double sided, 27 cm (60 and 100 grain)
312 230 500

Hard skin file stainless steel

23.5 cm; for sterile, individually packaged
polishing foils in coarse and medium*.
Especially important for callused skin tears -
e. g. in diabetics. Rapid change for each
new patient. The stainless steel handle is
easily sterilized.
Without polishing foils. Please order separately.
312 230 600

Polishing foils, 80 items.

Individually sterile packed in a
dispenser Grain size 60, coarse
312 230 601
Grain size 100, medium 312 230 602

Holding bench for callus file, prevents

scratches on the cabinet work area
312 230 700

Double hard skin file

with plastic handle 23 cm
SN 32 312 230 400



Secora hard skin plane, red



Right Side

KAPNOL ceramic hard skin plane 20 cm

312 100 900

Secora hard skin plane, red

steriliable, 17.3 cm
red, 0.50 mm 312 101 200
anthracite, 0.25 mm 312 101 300

KAPNOL hard skin lengthwise plane,

14.5 cm sterilisable
M 4 312 100 400

Lengthwise hard skin plane M 4



Special instruments, Files and Rasps

Rasps and files

Removal of hard skin and calluses.

Order No.

KAPNOL flat hard skin rasp

KAPNOL flat hard skin rasp

19 cm

312 100 500



Saphir nail file

Saphir nail file,

17.4 cm, with plastic handle, pointed

312 230 200



Sand paper file

Sand paper file,

17 cm

6 items

318 020 000



Nail file, notched

Nail file, notched, 15 cm

Stainless steel, pointed

HF 566

312 226 600



Nail polishing file

polishing file,

4.5 cm with 3 different grains:
For scratches and ridges on the nail

For matt shine on the nail plate
For shiny nail surfaces

312 250 000



Orthonyxy Requirements



Brace holder




Wire clipper



Orthonyxy three stage clipper



Loop bending clipper

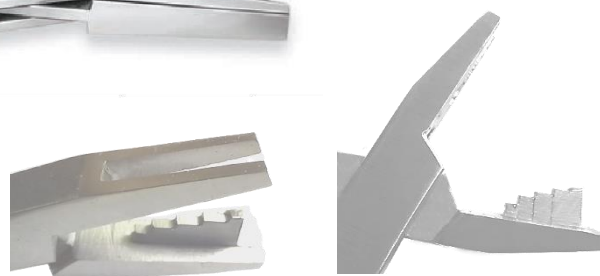



1 Round hammer



6 Wire cutter clipper R



7 Tuning fork

- All work processes
- Wire clipper DP O31R, 15 cm, Stainless steel, KAPNOL, essential
- 2 For creating nail braces and for shaping

and adjusting the brace

- Orthonyxy three stage clipper, 14 cm with blade, 3 differently sized loop bending devices and pointed mouth

- 4 Loop bending clipper, 14 cm, 7 stages

- 5 Spring steel wire
- Round hammer, 30 cm, for flattening the nail brace hook
- Wire cutter clipper R, 13 cm,

313 103 100

313 124 100

313 260 300

313 128 000

Spring steel wire

Sterilisable and stainless with hard metal insets on the blades

Rydel-Seiffer tuning fork, 23.5 cm, for measuring vibrations to establish the presence of polyneuropathy in diabetics for flattening nail braces

Spring steel wire 0.5 mm x 50 m for Ross Fraser braces



312 260 000

313 150 500


Brief application instructions



1. Measuring for the brace sections




2. Shortening the bent brace ends

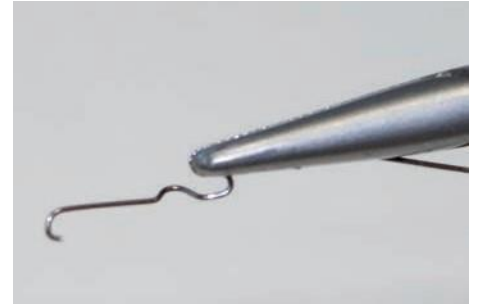

3. Shortened brace sections



4. Correctly applied brace sections

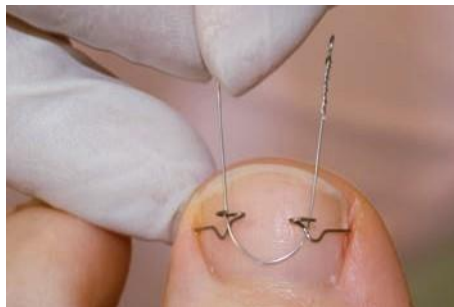

5. Section link to the GTO loop




6. Results after tightening the loop



Winding hook (313 270 000)



7. After removing all wire ends

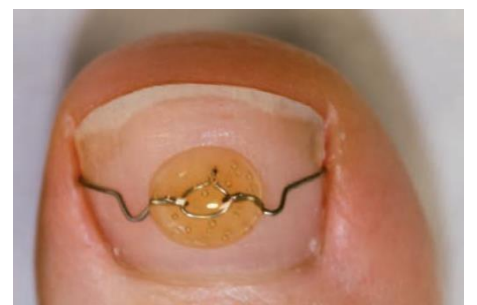

8. Final results

Wire thickness 0.3 mm (5 pcs.)
Wire thickness 0.4 mm (5 pcs.)

Order No. MD

313 340 100
313 340 200

Wire thickness 0.3 mm (5 pcs.)
Wire thickness 0.4 mm (5 pcs.)

313 340 300
313 340 400

Winding hook, stainless steel
Length 15.5 cm, special instrument for
applying the GTO nail brace.
The winding hook is used to twist the two
hooks on the brace together with the loop.

313 270 000

Using the VHO Osthold PERFECT brace



Two hooked steel wire springs are attached painlessly to the edges of the nail and pulled together with a loop. This pulls the nail edges out of the inflamed nail fold. There is no need for costly silicon impressions and plaster models, and no need for waiting.

1 The VHO brace consists of 3 parts: a, b and c. a and b are two tiny pre-formed levers on the right and left of a steel wire spring. In the middle is a shape designed to be guided by 3 fingers of the left hand. a and b are then firmly joined to each other on the nail plate by turning the loop, which is then cut off and glued in place.

2 The small lever is fitted precisely to one side of the nail. The wire must be adjusted depending on the strength and curvature of the nail and the hook h is shortened. It is fixed in place using high performance adhesive and cut off.




3 The other side is fitted in the same manner. This is not cut off.

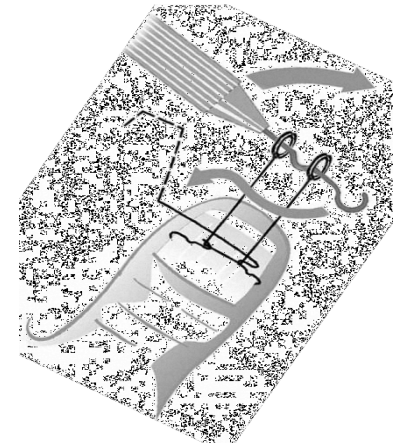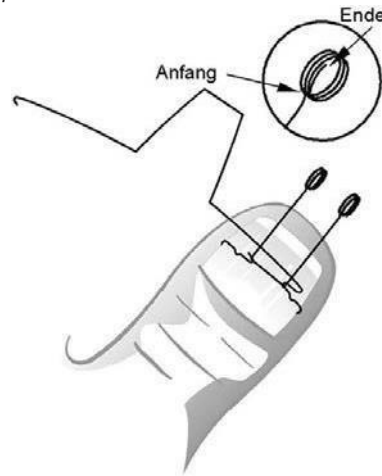

4 The loop is selected by visual judgment and wound round the small protrusions. The eye of the loop should be facing you. The loops are pressed together.

5 The winder is threaded through both eyes and positioned vertically.

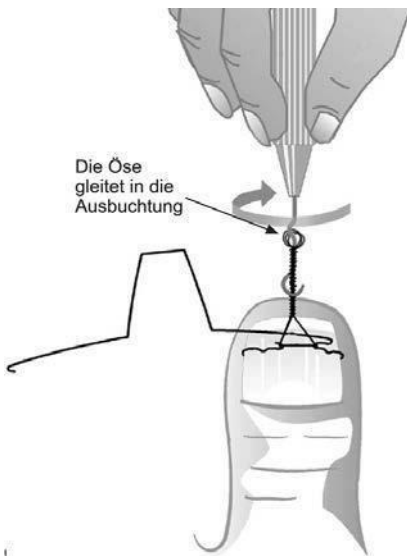

6 By turning the winder (clockwise) the three parts of the VHO bound together.



7 All wire parts that stick out are cut off as short as possible.

8 All the sharp points are covered by a thick drop of plastic nail material. If the nails are very curved, the wire parts can be covered with plastic nail after being placed on one side, either to the right or to the left





KAPNOL MEDICAL INSTRUMENTS



Winding hooks



2 Dolphin pliers



3 Round/pointed pliers

L

Winding hooks, stainless steel, rust-free
Length 15.5 cm, special instruments
for placing the VHO Osthold brace
The two small hooks on the brace are
wound together with the loop using the
winding hook.

Order No. MD
[313 270 000](#)

Stainless steel, **KAPNOL Dolphin Pliers**
Length 15.5 cm for bending the
necessary pulling power on the
tail and for adjusting the angle of the brace.

[313 260 600](#)

Round/pointed pliers, chrome
length 11.5 cm Custom production
adjust the hook and bring it into the
So called "submissive position".
wire cutter pliers.

[313 260 500](#)



DR. KAPNOL

PRC GOLF TOURNAMENT

HIGH QUALITY &
AFFORDABLE PRICE

KAPNOL MEDICAL INSTRUMENTS

Web: www.kapnol.com

Email: info@kapnol.com

Phone[#] (+92) 333 874 1182
(+92) 331 721 6007